

**Profenofos**

**Profenofos EC (A8591B) - *Salmonella Typhimurium* and  
*Escherichia Coli* Reverse Mutation Assay**

**Final Report**

**TEST GUIDELINE(S):** OECD 471 (2020)

**AUTHOR(S):** Dipl. Biol. Andrea Sokolowski

**COMPLETION DATE:** 18 November 2021

**PERFORMING LABORATORY:** ICCR-Roßdorf GmbH  
In den Leppsteinswiesen 19  
64380 Rossdorf, Germany

**LABORATORY PROJECT ID:** Report Number: 2180100  
Study Number: 2180100  
Task Number: TK0612496

**SPONSOR(S):** Syngenta Ltd.  
Jealott's Hill International Research Centre  
Bracknell, Berkshire RG42 6EY, United Kingdom

## STATEMENT OF DATA CONFIDENTIALITY CLAIMS

**The Following Statement Applies To The United States of America:**

### STATEMENT OF NO DATA CONFIDENTIALITY CLAIMS UNDER SPECIFIED FIFRA PROVISIONS

No claim of confidentiality, on any basis whatsoever, is made for any information contained in this document. I acknowledge that information not designated as within the scope of FIFRA sec. 10(d)(1)(A), (B), or (C) and which pertains to a registered or previously registered pesticide is not entitled to confidential treatment and may be released to the public, subject to the provisions regarding disclosure to multinational entities under FIFRA 10(g).

Company: Syngenta Crop Protection, LLC  
410 Swing Road  
Post Office Box 18300  
Greensboro, NC 27419-8300 USA

Submitter: \_\_\_\_\_ Date: \_\_\_\_\_

Syngenta is the owner of this information and data. Syngenta has submitted this material to the United States Environmental Protection Agency specifically under the provisions contained in FIFRA as amended and, hereby, consents to use and disclosure of this material by EPA according to FIFRA. In submitting this material to EPA according to method and format requirements contained in PR Notice 2011-3, we do not waive any protection or right involving this material that would have been claimed by the company if this material had not been submitted to the EPA, nor do we waive any protection or right provided under FIFRA Section 3 (concerning data exclusivity and data compensation) or FIFRA Section 10(g) (prohibiting disclosure to foreign and multinational pesticide companies or their agents).

## GOOD LABORATORY PRACTICE COMPLIANCE STATEMENT

This study performed in the test facility of ICCR-Roßdorf GmbH, In den Leppsteinswiesen 19, 64380 Rossdorf, Germany was conducted in compliance with Good Laboratory Practice Regulations:

Chemikaliengesetz (Chemicals Act) of the Federal Republic of Germany, "Anhang 1" (Annex 1), in its currently valid version


OECD Principles of Good Laboratory Practice, (as revised in 1997), ENV/MC/CHEM(98)17

EC Commission Directive 2004/10/EC

These procedures are compatible with Good Laboratory Practice regulations specified by regulatory authorities throughout the European Community, the United States (EPA and FDA), and Japan (MHW, MAFF, and METI), and other countries that are signatories to the OECD Mutual Acceptance of Data Agreement.

There were no circumstances that may have affected the quality or integrity of the study.

Dipl. Biol. Andrea Sokolowski  
Study Director Bacterial Systems

  
Date: 18 November 2021

Performing Laboratory:  
ICCR-Roßdorf GmbH  
In den Leppsteinswiesen 19  
64380 Rossdorf, Germany

To be completed for USA EPA submission only:  
Representative of Submitter/Sponsor:

_____	_____
	Date
Submitter/Sponsor: Syngenta Crop Protection, LLC 410 Swing Road Post Office Box 18300 Greensboro, NC 27419-8300 USA	

## **FLAGGING STATEMENT**

This page is intentionally left blank. It will be replaced by an appropriate Flagging statement by the Sponsor.

## QUALITY ASSURANCE STATEMENT

ICCR Study Number: 2180100  
Test substance: Profenofos EC (A8591B)  
Study director: Dipl. Biol. Andrea Sokolowski  
Study Title: Profenofos EC (A8591B) -  
*Salmonella Typhimurium* and  
*Escherichia Coli* Reverse Mutation Assay

Study based activities at the Test Facility ICCR-Roßdorf GmbH were audited and inspected.  
The details of these audits and inspections are given below.

Type of Inspection	Date(s) of Inspection	Date Reporting to Study Director, Test Facility Management
Study Plan Verification	06 August 2021	06 August 2021
Process – based Test system preparation and application	08 September 2021	08 September 2021
Report Audit	28 October 2021	29 October 2021

General facilities and activities where this study was conducted were inspected on an annual basis and results are reported to the relevant responsible person and Management.

The statement is to confirm, that this report reflects the raw data.



H. Pilawa

Quality Assurance Auditor  
ICCR-Roßdorf GmbH

18 November 2021

Date

## PROJECT STAFF SIGNATURE

Study Director

Dipl. Biol. Andrea Sokolowski

A handwritten signature consisting of the letters 'A' and 'S' in a cursive style, written over a horizontal dotted line.

Date: 18 November 2021

## GENERAL INFORMATION

### Contributors

The following contributed to this report in the capacities indicated:

Name	Title
Dipl. Biol. Andrea Sokolowski	Study Director
Dr. Markus Schulz	Test Facility Management
Frauke Hermann	Head of Quality Assurance Unit
Carolina Vaccari	Syngenta Study Manager

### Study Dates

Study initiation date:	09 August 2021
Experimental start date:	11 August 2021
Experimental completion date:	08 September 2021

### Deviations from the Guidelines

None

### Retention of Samples

None

### Performing Laboratory Test Substance Reference Number

S 2177311

### Other

ICCR-Roßdorf GmbH will archive:

Records and documentation relating to this study will be maintained in the archives of ICCR-Roßdorf GmbH for a period of 4 years from the date on which the Study Director signs the final report. This will include electronic and paper raw data, and report that support the reconstruction of the study.

At termination of the aforementioned period, the records and documentation will be transferred to the GLP compliant Archive of Rhenus Archiv Services GmbH, Frankfurt am Main for further archiving up to a total archiving period of 15 years.

A sample of the test substance will not be archived.

ICCR Roßdorf GmbH will retain in its archive a copy of the study plan and final report, and any amendments indefinitely.

### Deviations from the study plan

There were no deviations (unplanned changes) from the study plan.

**Distribution of the report**

Sponsor	2 × electronic copy (1 × pdf-file, 1 × Word-file)
Study Director	1 × (original)



## TABLE OF CONTENTS

<b>STATEMENT OF DATA CONFIDENTIALITY CLAIMS</b>	<b>2</b>
<b>GOOD LABORATORY PRACTICE COMPLIANCE STATEMENT</b>	<b>3</b>
<b>FLAGGING STATEMENT</b>	<b>4</b>
<b>QUALITY ASSURANCE STATEMENT</b>	<b>5</b>
<b>PROJECT STAFF SIGNATURE</b>	<b>6</b>
<b>GENERAL INFORMATION</b>	<b>7</b>
<b>TABLE OF CONTENTS</b>	<b>9</b>
<b>1.0 EXECUTIVE SUMMARY</b>	<b>11</b>
1.1 Study Design .....	11
1.2 Results .....	11
1.3 Conclusion.....	11
<b>2.0 INTRODUCTION</b>	<b>12</b>
2.1 Purpose.....	12
2.2 Test Guideline(s).....	12
<b>3.0 MATERIALS AND METHODS</b>	<b>13</b>
3.1 Test Substance.....	13
3.2 Controls .....	14
3.2.1 Negative controls .....	14
3.2.2 Positive control substances .....	14
3.3 Experimental Design.....	15
3.3.1 Characterisation of the <i>Salmonella typhimurium</i> and <i>E. coli</i> strains .....	15
3.3.2 Storage.....	16
3.3.3 Precultures.....	16
3.3.4 Selective agar .....	16
3.3.5 Overlay agar .....	16
3.4 Mammalian Microsomal Fraction S9 Mix.....	16
3.4.1 S9 (Preparation by ICCR-Roßdorf GmbH).....	17
3.4.2 S9 mix .....	17
3.5 Pre-Experiment for Cytotoxicity.....	17
3.6 Concentration Selection .....	18
3.7 Experimental Performance.....	18
3.8 Data Evaluation.....	19

3.8.1	Data recording .....	19
3.8.2	Acceptability of the assay .....	19
3.8.3	Evaluation of results.....	19
3.8.4	Biometry.....	20
3.9	Major Computerized System.....	20
<b>4.0</b>	<b>RESULTS AND DISCUSSION</b>	<b>21</b>
<b>5.0</b>	<b>CONCLUSIONS</b>	<b>22</b>
<b>6.0</b>	<b>REFERENCES</b>	<b>23</b>
	<b>TABLES SECTION</b>	<b>24</b>
TABLE 1	Summary of Results Pre-Experiment/Experiment I.....	25
TABLE 2	Summary of Results Experiment II.....	26
TABLE 3	Pre-Experiment/Experiment I: 2180100 VV Plate Incorporation Without Metabolic Activation.....	27
TABLE 4	Pre-Experiment/Experiment I: 2180100 VV Plate Incorporation With Metabolic Activation.....	29
TABLE 5	Experiment II: 2180100 HV2 Pre Incubation Without Metabolic Activation.....	31
TABLE 6	Experiment II: 2180100 HV2 Pre Incubation With Metabolic Activation.....	33
	<b>APPENDICES SECTION</b>	<b>35</b>
APPENDIX 1	Historical Control Data .....	36
APPENDIX 2	Copy of GLP Certificate .....	37
APPENDIX 3	Certificate of S9 .....	38
APPENDIX 4	Certificate of Analysis.....	39

## **1.0 EXECUTIVE SUMMARY**

### **1.1 Study Design**

This study was performed to investigate the potential of profenofos EC (A8591B) to induce gene mutations in the plate incorporation test (Experiment I) and the pre-incubation test (Experiment II) using the *Salmonella typhimurium* (*S. typhimurium*) strains TA1535, TA1537, TA98, and TA100, and the *Escherichia coli* (*E. coli*) strains WP2 *uvrA* (pKM101) and WP2 (pKM101).

### **1.2 Results**

The plates incubated with the test substance showed normal background growth up to the maximal concentration of 5000 µg/plate with and without S9 mix in all strains used.

Cytotoxic effects, evident as a reduction in the number of revertants (below the indication factor of 0.5), occurred in strains TA100 and WP2 (pKM101) in the presence and absence of S9 mix and in strain WP2 *uvrA* (pKM101) in the presence of S9 mix in both experiments.

No relevant increase in revertant colony numbers of any of the six tester strains was observed following treatment with profenofos EC (A8591B) at any concentration, neither in the presence nor absence of metabolic activation (S9 mix). There was also no observed tendency of higher mutation rates with increasing concentrations in the range below the generally acknowledged border of biological relevance.

Appropriate reference mutagens were used as positive controls, which showed a distinct increase of induced revertant colonies consistent with the laboratory's historical control data demonstrated the sensitivity of the test system and the efficacy of the S9 mix. Each batch of S9 was also tested with 2 pro-mutagens, benzo(a)pyrene and 2-aminoanthracene.

### **1.3 Conclusion**

In conclusion, it can be stated that during the described mutagenicity tests and under the experimental conditions reported, profenofos EC (A8591B) did not induce gene mutations by base pair changes or frameshifts in the genome of the strains used.

Therefore, profenofos EC (A8591B) is considered to be negative (i.e. non-mutagenic) in the *Salmonella typhimurium* and *Escherichia coli* reverse mutation assay.

## 2.0 INTRODUCTION

### 2.1 Purpose

These experiments were performed to assess the potential of the test substance to induce gene mutations by means of the *S. typhimurium* and *E. coli* reverse mutation assay. Experiment I was performed as a plate incorporation assay. Since a negative result was obtained in this experiment, Experiment II was performed as a pre-incubation assay.

The most widely used assays for detecting gene mutations are those using bacteria (1). They are relatively simple and rapid to perform, and give reliable data on the ability of an agent to interact with DNA and produce mutations.

Reverse mutation assays determine the frequency with which an agent reverses or suppresses the effect of the forward mutation. The genetic target presented to an agent is therefore small, specific and selective. Several bacterial strains, or a single strain with multiple markers are necessary to assure reliable detection of mutagens that may be specific to one tester strain or locus. The reversion of bacteria from growth-dependence on a particular amino acid to growth in the absence of that amino acid (reversion from auxotrophy to prototrophy) is the most widely used marker.

The *S. typhimurium* histidine (his) and the *E. coli* tryptophan (trp) reversion system measures  $\text{his}^- \rightarrow \text{his}^+$  and  $\text{trp}^- \rightarrow \text{trp}^+$  reversions, respectively. The *S. typhimurium* and *E. coli* strains are constructed to differentiate between base pair (TA1535, TA100, WP2 *uvrA* (pKM101), and WP2 (pKM101)) and frameshift (TA1537, TA98) mutations.

According to the direct plate incorporation and pre-incubation method the bacteria are exposed to the test substance with and without metabolic activation and plated on selective medium. After a suitable period of incubation, revertant colonies are counted.

To establish a concentration response effect at least seven concentrations with adequately spaced intervals were tested. The maximum concentration was 5000 µg/plate.

To validate the test, reference mutagens were tested in parallel to the test substance.

### 2.2 Test Guideline(s)

This study followed the procedures indicated by the following internationally accepted guideline and recommendations:

“Ninth Addendum to OECD Guidelines for Testing of Chemicals”, Section 4, No. 471:  
“Bacterial Reverse Mutation Test”, corrected June 26, 2020

### 3.0 MATERIALS AND METHODS

#### 3.1 Test Substance

Information as provided by the Sponsor.

Identification:	Profenofos EC (A8591B)
Batch:	RAN001-099-019
Content of Profenofos:	73.12 % w/w corresponding to 969.03 g/L
Appearance:	Liquid, colourless to yellowish*
Recertification Date:	06 May 2023
Storage Conditions:	At room temperature
Stability in Solvent:	Not indicated by the Sponsor

\*determined by ICCR-Roßdorf staff

The test substance concentrations were not adjusted for the content of profenofos.

On the day of the experiment (immediately before use), the test substance was dissolved in dimethylsulfoxide (DMSO, purity > 99%). The solvent was chosen as the best suitable solvent compared to water and ethanol, according to its solubilisation properties and its relative non-toxicity to the bacteria (2).

All formulations were prepared freshly before treatment and used within two hours of preparation. The formulation was assumed to be stable for this period unless specified otherwise by the Sponsor.

## 3.2 Controls

### 3.2.1 Negative controls

Concurrent untreated and solvent controls were performed.

### 3.2.2 Positive control substances

#### Without metabolic activation

Strains:	TA1535, TA100
Name:	Sodium azide, (NaN <sub>3</sub> )
Supplier:	Sigma-Aldrich, 82024 Taufkirchen, Germany
Batch No.:	STBJ7813
Purity:	≥ 99%
Dissolved in:	Deionised water
Concentration:	10 µg/plate
Strains:	TA1537, TA98
Name:	4-nitro-o-phenylene-diamine, (4-NOPD)
Supplier:	Sigma-Aldrich, 82024 Taufkirchen, Germany
Batch No.:	MKBM 5257V
Purity:	≥ 98%
Dissolved in:	DMSO (purity >99 %, Fisher Leics LE11 5RG, United Kingdom)
Concentration:	10 µg/plate in strain TA 98, 50 µg/plate in strain TA 1537
Strains:	WP2 <i>uvrA</i> (pKM101), WP2 (pKM101)
Name:	Methyl methane sulfonate, (MMS)
Supplier:	Sigma-Aldrich, 82024 Taufkirchen, Germany
Batch No.:	MKCL 6261
Purity:	≥ 99%
Dissolved in:	Deionised water
Concentration:	2.0 µL/plate

#### With metabolic activation

Strains:	TA1535, TA1537, TA98, TA100, WP2 <i>uvrA</i> (pKM101), WP2 (pKM 101)
Name:	2-aminoanthracene, (2-AA)
Supplier:	Sigma-Aldrich, 82024 Taufkirchen, Germany
Batch No.:	STBG 0630V
Purity:	≥ 96%
Dissolved in:	DMSO (purity > 99 %, Fisher Leics LE11 5RG, United Kingdom)
Concentration:	2.5 µg/plate (TA1535, TA1537, TA98, TA100), 10 µg/plate (WP2 <i>uvrA</i> (pKM101), WP2 (pKM101))

The stability of the positive control substances in solution is unknown but a mutagenic response in the expected range is sufficient evidence of biological activity.

### 3.3 Experimental Design

#### 3.3.1 Characterisation of the *Salmonella typhimurium* and *E. coli* strains

The histidine dependent strains are derived from *S. typhimurium* strain LT2 through mutations in the histidine locus. Additionally, due to the "deep rough" (*rfa*<sup>-</sup>) mutation they possess a faulty lipopolysaccharide envelope which enables substances to penetrate the cell wall more easily. A further mutation causes a reduction in the activity of an excision repair system. The last alteration includes mutational processes in the nitrate reductase and biotin genes produced in a UV-sensitive area of the gene named *uvrB*<sup>-</sup>. In the strains TA98 and TA100 the R-factor plasmid pKM101 carries the ampicillin resistance marker (3).

Strain WP2 (4) and its derivatives all carry the same defect in one of the genes for tryptophan biosynthesis. Tryptophan-independent (*Trp*<sup>+</sup>) mutants (revertants) can arise either by a base change at the site of the original alteration or by a base change elsewhere in the chromosome so that the original defect is suppressed. This second possibility can occur in several different ways so that the system seems capable of detecting all types of mutagen which substitute one base for another. Additionally, the *uvrA* derivative is deficient in the DNA repair process (excisable repair damage). Such a repair-deficient strain may be more readily mutated by agents. The *E. coli* strains WP2 *uvrA* (pKM101) and WP2 (pKM101) are constructed by introduction of the R-factor plasmid pKM101.

When summarized, the mutations of the *S. typhimurium* and *E. coli* strains used in this study can be described as follows:

Strains	Genotype	Type of mutations indicated
<b><i>Salmonella typhimurium</i></b>		
TA1537	<i>his</i> C 3076; <i>rfa</i> <sup>-</sup> ; <i>uvrB</i> <sup>-</sup>	frame shift mutations
TA98	<i>his</i> D 3052; <i>rfa</i> <sup>-</sup> ; <i>uvrB</i> <sup>-</sup> ; R-factor	" "
TA1535	<i>his</i> G 46; <i>rfa</i> <sup>-</sup> ; <i>uvrB</i> <sup>-</sup>	base-pair substitutions
TA100	<i>his</i> G 46; <i>rfa</i> <sup>-</sup> ; <i>uvrB</i> <sup>-</sup> ; R-factor	" "
<b><i>Escherichia coli</i></b>		
WP2 <i>uvrA</i> (pKM101)	<i>trp</i> E 56 <i>uvrA</i> <sup>-</sup> ; R-factor	base-pair substitutions and others
WP2 (pKM101)	<i>trp</i> E 56; R-factor	" "

Regular checking of the properties of the *S. typhimurium* and *E. coli* strains regarding the membrane permeability and ampicillin resistance; UV sensitivity, and amino acid requirement as well as normal spontaneous mutation rates is performed by ICCR-Roßdorf GmbH according to Ames *et al.* (5), Maron and Ames (3), and Mortelmans and Riccio (7). In this way it is ensured that the experimental conditions set down by Ames are fulfilled.

The bacterial strains TA1535, TA1537, TA98, TA100, WP2 *uvrA* (pKM101), and WP2 (pKM101) were obtained from Trinova Biochem GmbH (35394 Gießen, Germany).

### 3.3.2 Storage

The strain cultures were stored as stock cultures in ampoules with nutrient broth + 5 % DMSO (Fisher Leics, LE11 5RG, United Kingdom) in liquid nitrogen.

### 3.3.3 Precultures

The thawed bacterial suspension was transferred into 250 mL Erlenmeyer flasks containing nutrient medium (50 mL). A solution of ampicillin (50 µL, 25 µg/mL) was added to the strains TA98, TA100, WP2 *uvrA* (pKM101), and WP2 (pKM101). This nutrient medium contains per liter:

8 g Nutrient Broth (MERCK, 64293 Darmstadt, Germany)  
5 g NaCl (MERCK, 64293 Darmstadt, Germany)

The bacterial cultures were incubated in a shaking water bath for 4 hours at 37° C. The optical density of the bacteria was determined by absorption measurement and the obtained values indicated that the bacteria were harvested at the late exponential or early stationary phase ( $10^8$ - $10^9$  cells/mL).

### 3.3.4 Selective agar

Plates with selective agar (without Histidine/Tryptophan) were used.

### 3.3.5 Overlay agar

The overlay agar contained per litre:

for *Salmonella* strains:

7.0 g Agar Agar\*

6.0 g NaCl\*

10.5 mg L-Histidine×HCl×H<sub>2</sub>O\*

12.2 mg Biotin\*

for *Escherichia coli* strains:

7.0 g Agar Agar\*

6.0 g NaCl\*

10.2 mg Tryptophan\*

\* (MERCK, 64293 Darmstadt, Germany)

Sterilisations were performed at 121° C in an autoclave.

## 3.4 Mammalian Microsomal Fraction S9 Mix

The bacteria used in this assay do not possess the enzyme systems which, in mammals, are known to convert promutagens into active DNA damaging metabolites. In order to overcome this major drawback an exogenous metabolic system is added in the form of mammalian microsome enzyme activation mixture.



### 3.4.1 S9 (Preparation by ICCR-Roßdorf GmbH)

Phenobarbital/ $\beta$ -naphthoflavone induced rat liver S9 was used as the metabolic activation system. The S9 was prepared from male Wistar rats (RjHan:WI; weight approx. 220 – 320 g, Janvier Labs, 53941 Saint-Berthevin Cedex, France) induced by peroral administration of 80 mg/kg b.w. phenobarbital (Sigma-Aldrich Chemie GmbH, 82024 Taufkirchen, Germany) and by peroral administrations of  $\beta$ -naphthoflavone (Acros Organics, 2440 Geel, Belgium) each, on three consecutive days. The livers were prepared 24 hours after the last treatment. The S9 fractions were produced by dilution of the liver homogenate with a KCl solution (1+3 parts) followed by centrifugation at 9000 g. Aliquots of the supernatant were frozen and stored in ampoules at  $-80^{\circ}\text{C}$ . Small numbers of the ampoules can be kept at  $-20^{\circ}\text{C}$  for up to one week. Each batch of S9 mix is routinely tested with 2-aminoanthracene as well as benzo[a]pyrene (Appendix 3).

The protein concentration in the S9 preparation was 30.2 mg/mL (lot no. 200521) in both experiments.

### 3.4.2 S9 mix

Before the experiment an appropriate quantity of S9 supernatant was thawed and mixed with S9 cofactor solution. The amount of S9 supernatant was 10% v/v in the S9 mix. Cofactors were added to the S9 mix to reach the following concentrations in the S9 mix:

8 mM  $\text{MgCl}_2$   
33 mM KCl  
5 mM Glucose-6-phosphate  
4 mM NADP

in 100 mM sodium-ortho-phosphate-buffer, pH 7.4.

During the experiment the S9 mix was stored in an ice bath. The S9 mix preparation was performed according to Ames *et al.*(5).

## 3.5 Pre-Experiment for Cytotoxicity

To evaluate the cytotoxicity of the test substance a pre-experiment was performed with all strains. Eight concentrations were tested for cytotoxicity and mutation induction each with three replicate plates. The experimental conditions in this pre-experiment are described in section 3.7 (plate incorporation test).

Cytotoxicity of the test substance results in a reduction in the number of spontaneous revertants (below a factor of 0.5) or a clearing of the bacterial background lawn.

The pre-experiment is reported as the Main Experiment I since the criteria mentioned in Section 3.8.2 Acceptability of the Assay were met.

### 3.6 Concentration Selection

In the pre-experiment the concentration range of the test substance was 3 - 5000 µg/plate. The pre-experiment is reported as Experiment I. Since minor cytotoxic effects and precipitation of the test substance occurred in Experiment I, seven concentrations were tested in Experiment II. 5000 µg/plate was chosen as the maximal concentration in Experiment II.

The concentration range included two logarithmic decades. The following concentrations were tested in Experiment II:

10; 33; 100; 333; 1000; 2500; and 5000 µg/plate

### 3.7 Experimental Performance

For each strain and concentration including the controls, three plates were used.

The following materials were mixed in a test tube and poured onto the selective agar plates:

- 100 µL Test solution at each concentration, solvent (negative control) or reference mutagen solution (positive control),
- 500 µL S9 mix (for test with metabolic activation) or S9 mix substitution buffer\* (for test without metabolic activation),
- 100 µL Bacteria suspension (cf. test system, pre-culture of the strains; OD = 1.0 - 1.2; wavelength = 500 nm; approx.  $8 \times 10^8$  cells/mL),
- 2000 µL Overlay agar

For the pre-incubation method test solution (100 µL) (solvent or reference mutagen solution (positive control)), S9 mix / S9 mix substitution buffer\* (500 µL) and bacteria suspension (100 µL) were mixed in a test tube and incubated at  $37^\circ \text{C} \pm 1.5^\circ \text{C}$  for 60 minutes. After pre-incubation overlay agar (2.0 mL,  $45^\circ \text{C}$ ) was added to each tube. The mixture was poured on selective agar plates.

After solidification the plates were incubated upside down for 72 hours at  $37^\circ \text{C} \pm 1.5^\circ \text{C}$  in the dark, plates were then stored at  $4^\circ \text{C}$  until counted (6).

In parallel to each test a sterile control of the test substance was performed and documented in the raw data. Therefore, stock solution (100 µL) and S9 mix / S9 mix substitution buffer\* (500 µL) were mixed with overlay agar (2.0 mL) and poured on minimal agar plates.

\* Substitution buffer: 7 parts of the 100 mM sodium-ortho-phosphate-buffer pH 7.4 with 3 parts of KCl solution 0.15 M

## **3.8 Data Evaluation**

### **3.8.1 Data recording**

The colonies were counted using a Petri Viewer with the software program Ames Study Manager (see section 3.9, Major computerized systems). The evaluation unit was connected to a PC with printer to print out the individual values, the means from the plates for each concentration together with standard deviations and enhancement factors as compared to the spontaneous reversion rates (see tables of results). The print outs are kept with the raw data. Due to precipitation of the test item some test groups were scored manually (as indicated on data tables).

### **3.8.2 Acceptability of the assay**

The *Salmonella typhimurium* and *Escherichia coli* reverse mutation assay is considered acceptable if it meets the following criteria:

- regular background growth in the negative and solvent control
- the spontaneous reversion rates in the negative and solvent control are in the range of the historical data
- the positive control substances should produce an increase above the threshold of twofold (strains TA 98, TA 100, WP2 uvrA (pKM101, and WP2 (pKM101))) or threefold (strains TA 1535 and TA 1537) the revertant colony count of the corresponding solvent control;
- a minimum of five analysable concentrations should be present with at least four showing no signs of toxic effects, evident as a reduction in the number of revertants below the indication factor of 0.5.

### **3.8.3 Evaluation of results**

A test substance is considered as a mutagen if a biologically relevant increase in the number of revertants of twofold or above (strains TA 98, TA 100, WP2 uvrA (pKM101), and WP2 (pKM101)) or of threefold or above (strains TA 1535 and TA 1537) the spontaneous mutation rate of the corresponding solvent control is observed (1).

A concentration dependent increase is considered biologically relevant if the threshold is reached or exceeded at more than one concentration (6).

An increase of revertant colonies equal or above the threshold at only one concentration is judged as biologically relevant if reproduced in an independent second experiment.

A concentration dependent increase in the number of revertant colonies below the threshold is regarded as an indication of a mutagenic potential if reproduced in an independent second experiment. However, whenever the colony counts remain within the historical range of negative and solvent controls, such an increase is not considered biologically relevant.

#### **3.8.4 Biometry**

According to the OECD guideline 471, a statistical analysis of the data is not mandatory.

### **3.9 Major Computerized System**

Petri Viewer Sorcerer Colony Counter 3.0 (Instem, Suffolk IP33 3TA, UK) with the software program Ames Study Manager (v1.24) and Ames Archive Manager (v1.01).

## 4.0 RESULTS AND DISCUSSION

The test substance, profenofos EC (A8591B), was assessed for its potential to induce gene mutations in the plate incorporation test (Experiment I) and the pre-incubation test (Experiment II) using *S. typhimurium* strains TA1535, TA1537, TA98, and TA100, and the *E. coli* strains WP2 (pKM101) and WP2 *uvrA* (pKM101).

In the pre-experiment the concentration range of the test substance was 3 - 5000 µg/plate. The pre-experiment is reported as Experiment I. Since minor cytotoxic effects and precipitation of the test substance occurred in Experiment I, seven concentrations were tested in Experiment II. 5000 µg/plate was chosen as the maximal concentration in Experiment II. This is the maximum concentration recommended in the OECD test guideline.

The assay was performed with and without liver microsomal activation. Each concentration, including the controls, was tested in triplicate. The concentration range included two logarithmic decades. The test substance was tested at the following concentrations:

Pre-Experiment/Experiment I: 3; 10; 33; 100; 333; 1000; 2500; and 5000 µg/plate

Experiment II: 10; 33; 100; 333; 1000; 2500; and 5000 µg/plate

The test substance precipitated in the overlay agar in the test tubes from 1000 to 5000 µg/plate. Precipitation of the test item in the overlay agar on the incubated agar plates was observed from 2500 to 5000 µg/plate.

The plates incubated with the test substance showed normal background growth up to the maximal dose of 5000 µg/plate with and without S9 mix in all strains used.

Cytotoxic effects, evident as a reduction in the number of revertants (below the induction factor of 0.5), were observed at the following concentrations (µg/plate):

Strain	Experiment I		Experiment II	
	without S9 mix	with S9 mix	without S9 mix	with S9 mix
TA1535	/	/	/	/
TA1537	/	/	/	/
TA98	/	/	/	/
TA100	2500 – 5000	2500 – 5000	2500 – 5000	2500 – 5000
WP2 (pKM101)	5000	5000	2500 – 5000	5000
WP2 <i>uvrA</i> (pKM101)	/	5000	/	5000

/ = no cytotoxic effects, evident as a reduction in the number of revertants (below the induction factor of 0.5)

No substantial increase in revertant colony numbers in any of the six tester strains was observed following treatment with profenofos EC (A8591B) at any concentration, neither in the presence nor absence of metabolic activation (S9 mix). There was also no tendency of higher mutation rates with increasing concentrations in the range below the generally acknowledged border of biological relevance.

Appropriate reference mutagens were used as positive controls. They showed a distinct increase in induced revertant colonies.

## **5.0 CONCLUSIONS**

In conclusion, it can be stated that during the described mutagenicity tests and under the experimental conditions reported, profenofos EC (A8591B) did not induce gene mutations by base pair changes or frameshifts in the genome of the strains used.

Therefore, profenofos EC (A8591B) is considered to be non-mutagenic in the *Salmonella typhimurium* and *Escherichia coli* reverse mutation assay.

## 6.0 REFERENCES

1. Hollstein, M., J. McCann, F.A. Angelosanto, and W.W. Nichols (1979)  
Short-term tests for carcinogens and mutagens  
Mutation Res. 65, 133-226
2. Maron, D.M., J. Katzenellenbogen, and B.N. Ames (1981)  
Compatibility of organic solvents with the Salmonella/Microsome Test  
Mutation Res. 88, 343-350
3. Maron, D.M. and B.N. Ames (1983)  
Revised methods for the Salmonella mutagenicity test  
Mutation Res. 113, 173-215
4. Green, M.H.L. and W.J. Muriel (1976)  
Mutagen Testing Using TRP<sup>+</sup> Reversion in Escherichia Coli  
Mutation Res. 38, 3-32
5. Ames, B.N., J. McCann, and E. Yamasaki (1977)  
Methods for detecting carcinogens and mutagens with the Salmonella/mammalian  
microsome mutagenicity test  
In: B.J. Kilbey et al. (Eds.) "Handbook of Mutagenicity Test Procedures" Elsevier,  
Amsterdam, 1-17
6. de Serres, F.J. and M.D. Shelby (1979)  
Recommendations on data production and analysis using the Salmonella/microsome  
mutagenicity assay  
Mutation Res. 64, 159-165
7. Mortelmans, K. and E.S. Riccio (2000)  
The bacterial tryptophan reverse mutation assay with Escherichia coli WP2  
Mutation Res. 455, 61-69

## **TABLES SECTION**



**TABLE 1                      Summary of Results Pre-Experiment/Experiment I**

Study Name: 2180100  
Experiment: 2180100 VV Plate  
Assay Conditions:

Study Code: ICCR 2180100  
Date Plated: 11.08.2021  
Date Counted: 18.08.2021

Metabolic Activation	Test Group	Concentration (per plate)	Revertant Colony Counts (Mean ±SD)					
			TA 1535	TA 1537	TA 98	TA 100	WP2 pKM101	WP2 uvrA pKM101
Without Activation	DMSO		13 ± 4	14 ± 1	22 ± 2	103 ± 17	296 ± 15	291 ± 13
	Untreated		13 ± 3	10 ± 3	34 ± 7	102 ± 9	328 ± 21	322 ± 16
	Profenofos EC (A8591B)	3 µg	17 ± 1	12 ± 3	25 ± 2	101 ± 6	287 ± 18	307 ± 7
		10 µg	16 ± 1	13 ± 2	24 ± 7	105 ± 10	317 ± 31	305 ± 24
		33 µg	14 ± 4	12 ± 3	22 ± 10	90 ± 6	282 ± 9	298 ± 9
		100 µg	13 ± 4	11 ± 2	23 ± 6	89 ± 6	272 ± 28	305 ± 11
		333 µg	12 ± 4	11 ± 3	23 ± 8	75 ± 14	271 ± 19	272 ± 12
		1000 µg	10 ± 1	8 ± 2	24 ± 4	66 ± 4	213 ± 26	232 ± 21
		2500 µg	11 ± 2 <sup>P</sup>	8 ± 1 <sup>P M</sup>	17 ± 2 <sup>P</sup>	44 ± 2 <sup>P</sup>	178 ± 16 <sup>P</sup>	200 ± 6 <sup>P</sup>
		5000 µg	10 ± 2 <sup>P M</sup>	8 ± 1 <sup>P M</sup>	17 ± 2 <sup>P M</sup>	23 ± 1 <sup>P M</sup>	112 ± 7 <sup>P</sup>	181 ± 8 <sup>P M</sup>
	NaN3	10 µg	1040 ± 74			1562 ± 96		
	4-NOPD	10 µg			721 ± 70			
	4-NOPD	50 µg		78 ± 4				
	MMS	2.0 µL					2749 ± 102	2964 ± 226
With Activation	DMSO		11 ± 3	19 ± 2	43 ± 13	101 ± 20	310 ± 34	366 ± 3
	Untreated Profenofos EC (A8591B)		12 ± 3	18 ± 3	49 ± 3	117 ± 13	340 ± 6	361 ± 25
		3 µg	14 ± 3	18 ± 3	53 ± 9	110 ± 9	326 ± 19	363 ± 7
		10 µg	14 ± 4	16 ± 5	53 ± 4	109 ± 11	300 ± 6	363 ± 9
		33 µg	13 ± 4	17 ± 3	57 ± 13	120 ± 3	279 ± 25	343 ± 21
		100 µg	13 ± 1	17 ± 2	50 ± 2	109 ± 6	297 ± 10	309 ± 30
		333 µg	12 ± 3	12 ± 3	33 ± 6	97 ± 6	254 ± 28	339 ± 25
		1000 µg	11 ± 4	11 ± 2	33 ± 11	76 ± 13	187 ± 15	276 ± 36
		2500 µg	9 ± 3 <sup>P</sup>	11 ± 2 <sup>P M</sup>	25 ± 8 <sup>P</sup>	40 ± 10 <sup>P</sup>	160 ± 21 <sup>P</sup>	235 ± 8 <sup>P</sup>
		5000 µg	10 ± 1 <sup>P M</sup>	11 ± 3 <sup>P M</sup>	26 ± 3 <sup>P M</sup>	29 ± 7 <sup>P M</sup>	112 ± 5 <sup>P M</sup>	163 ± 13 <sup>P M</sup>
	2-AA	2.5 µg	297 ± 8	385 ± 48	2350 ± 78	3440 ± 40		
	2-AA	10.0 µg					952 ± 70	1670 ± 45
Key to Positive Controls			Key to Plate Postfix Codes					
NaN3	sodium azide		P	Precipitate				
2-AA	2-aminoanthracene		M	Manual count				
4-NOPD	4-nitro-o-phenylene-diamine							
MMS	methyl methane sulfonate							

**TABLE 2                      Summary of Results Experiment II**

Study Name: 2180100  
Experiment: 2180100 HV2 Pre  
Assay Conditions:

Study Code: ICCR 2180100  
Date Plated: 01.09.2021  
Date Counted: 08.09.2021

Metabolic <u>Activation</u>	Test <u>Group</u>	Concen- tration (per plate)	Revertant Colony Counts (Mean ±SD)					
			<u>TA 1535</u>	<u>TA 1537</u>	<u>TA 98</u>	<u>TA 100</u>	<u>WP2</u> <u>pKM101</u>	<u>WP2 <i>uvrA</i></u> <u>pKM101</u>
Without Activation	DMSO		14 ± 2	12 ± 3	23 ± 2	105 ± 6	263 ± 31	275 ± 6
	Untreated		16 ± 5	10 ± 2	33 ± 6	94 ± 12	324 ± 25	314 ± 11
	Profenofos EC (A8591B)	10 µg	13 ± 3	10 ± 4	26 ± 5	103 ± 9	260 ± 5	268 ± 21
		33 µg	11 ± 1	10 ± 2	28 ± 9	93 ± 7	264 ± 11	258 ± 8
		100 µg	14 ± 5	12 ± 3	21 ± 6	82 ± 11	178 ± 23	250 ± 3
		333 µg	12 ± 3	11 ± 3	20 ± 5	84 ± 9	155 ± 10	256 ± 14
		1000 µg	8 ± 2	8 ± 3	23 ± 2	76 ± 8	119 ± 22	181 ± 12
		2500 µg	9 ± 2 <sup>P</sup>	8 ± 1 <sup>P</sup>	20 ± 0 <sup>P</sup>	38 ± 9 <sup>P</sup>	96 ± 18 <sup>P</sup>	179 ± 15 <sup>P</sup>
		5000 µg	9 ± 0 <sup>PM</sup>	6 ± 2 <sup>PM</sup>	13 ± 3 <sup>PM</sup>	18 ± 3 <sup>PM</sup>	64 ± 8 <sup>P</sup>	133 ± 3 <sup>P</sup>
	NaN3	10 µg	938 ± 44			1425 ± 143		
	4-NOPD	10 µg			777 ± 41			
	4-NOPD	50 µg		100 ± 14				
	MMS	2.0 µL					2538 ± 35	2205 ± 131
With Activation	DMSO		11 ± 2	16 ± 1	45 ± 10	101 ± 10	312 ± 27	373 ± 16
	Untreated		11 ± 1	13 ± 3	39 ± 10	99 ± 17	319 ± 22	392 ± 8
	Profenofos EC (A8591B)	10 µg	14 ± 2	14 ± 3	44 ± 13	104 ± 16	291 ± 14	369 ± 27
		33 µg	14 ± 3	13 ± 3	46 ± 4	114 ± 15	272 ± 30	360 ± 14
		100 µg	11 ± 1	13 ± 2	48 ± 5	103 ± 17	284 ± 30	356 ± 22
		333 µg	12 ± 3	12 ± 2	35 ± 7	86 ± 6	219 ± 27	358 ± 26
		1000 µg	14 ± 3	14 ± 2	28 ± 5	87 ± 9	170 ± 10	270 ± 18
		2500 µg	14 ± 3 <sup>P</sup>	13 ± 2 <sup>P</sup>	28 ± 2 <sup>P</sup>	39 ± 7 <sup>P</sup>	175 ± 28 <sup>P</sup>	221 ± 8 <sup>P</sup>
		5000 µg	9 ± 3 <sup>PM</sup>	10 ± 1 <sup>PM</sup>	21 ± 6 <sup>PM</sup>	25 ± 4 <sup>PM</sup>	119 ± 8 <sup>P</sup>	164 ± 6 <sup>PM</sup>
	2-AA	2.5 µg	313 ± 23	307 ± 43	2288 ± 268	2502 ± 175		
	2-AA	10.0 µg					959 ± 34	1452 ± 90
Key to Positive Controls			Key to Plate Postfix Codes					
NaN3	sodium azide			P	Precipitate			
2-AA	2-aminoanthracene			M	Manual count			
4-NOPD	4-nitro-o-phenylene-diamine							
MMS	methyl methane sulfonate							

**TABLE 3                      Pre-Experiment/Experiment I: 2180100 VV Plate Incorporation Without Metabolic Activation**

Study Name: 2180100  
Experiment: 2180100 VV Plate  
Assay Conditions:

Study Code: ICCR 2180100  
Date Plated: 11.08.2021  
Date Counted: 18.08.2021

Without metabolic activation						
Strain	Compound	Concentration per plate	Mean revertants per plate	Standard Deviation	Ratio treated / solvent	Individual revertant colony counts
TA 1535	Profenofos EC (A8591B)	3 µg	16.7	0.6	1.3	17, 16, 17
		10 µg	15.7	0.6	1.2	16, 16, 15
		33 µg	13.7	4.0	1.1	16, 16, 9
		100 µg	13.0	3.6	1.0	9, 16, 14
		333 µg	12.3	3.5	0.9	12, 9, 16
		1000 µg	9.7	1.2	0.7	11, 9, 9
		2500 µg	10.7	1.5	0.8	11 P, 12 P, 9 P
		5000 µg	9.7	2.1	0.7	8 P M, 12 P M, 9 P M
	DMSO		13.0	3.6		12, 10, 17
	Untreated		13.0	2.6		15, 14, 10
TA 1537	Profenofos EC (A8591B)	3 µg	12.3	2.5	0.9	12, 15, 10
		10 µg	13.3	2.1	0.9	11, 15, 14
		33 µg	11.7	2.9	0.8	10, 15, 10
		100 µg	10.7	1.5	0.7	11, 9, 12
		333 µg	11.0	2.6	0.8	14, 10, 9
		1000 µg	8.0	1.7	0.6	9, 9, 6
		2500 µg	8.3	1.2	0.6	7 P M, 9 P M, 9 P M
		5000 µg	7.7	0.6	0.5	8 P M, 7 P M, 8 P M
	DMSO		14.3	0.6		15, 14, 14
	Untreated		10.0	3.5		6, 12, 12
TA 98	Profenofos EC (A8591B)	3 µg	25.0	2.0	1.2	23, 27, 25
		10 µg	23.7	7.2	1.1	20, 19, 32
		33 µg	22.0	9.5	1.0	16, 33, 17
		100 µg	22.7	6.4	1.0	30, 19, 19
		333 µg	22.7	7.5	1.0	27, 27, 14
		1000 µg	23.7	4.0	1.1	23, 20, 28
		2500 µg	17.0	1.7	0.8	16 P, 19 P, 16 P
		5000 µg	16.7	2.1	0.8	15 P M, 19 P M, 16 P M
	DMSO		21.7	1.5		22, 23, 20
	Untreated		34.3	6.7		27, 40, 36

Key to Plate Postfix Codes

P      Precipitate  
M      Manual count

Study Name: 2180100  
 Experiment: 2180100 VV Plate  
 Assay Conditions:

Study Code: ICCR 2180100  
 Date Plated: 11.08.2021  
 Date Counted: 18.08.2021

**Without metabolic activation**

Strain	Compound	Concentration per plate	Mean revertants per plate	Standard Deviation	Ratio treated / solvent	Individual revertant colony counts
<b>TA 100</b>	<b>Profenofos EC (A8591B)</b>	3 µg	101.3	5.5	1.0	105, 95, 104
		10 µg	104.7	9.9	1.0	98, 116, 100
		33 µg	89.7	5.7	0.9	96, 88, 85
		100 µg	89.0	6.0	0.9	83, 89, 95
		333 µg	75.0	14.0	0.7	65, 69, 91
		1000 µg	65.7	4.2	0.6	61, 69, 67
		2500 µg	44.3	1.5	0.4	46 P, 43 P, 44 P
		5000 µg	23.3	0.6	0.2	23 P M, 24 P M, 23 P M
	<b>DMSO</b>		102.7	17.2		115, 83, 110
	<b>Untreated</b>		102.3	9.1		101, 112, 94
<b>WP2 pKM101</b>	<b>Profenofos EC (A8591B)</b>	3 µg	287.3	17.6	1.0	307, 282, 273
		10 µg	316.7	31.1	1.1	282, 326, 342
		33 µg	281.7	8.5	1.0	272, 288, 285
		100 µg	271.7	27.6	0.9	291, 284, 240
		333 µg	271.3	19.0	0.9	253, 291, 270
		1000 µg	212.7	25.6	0.7	215, 237, 186
		2500 µg	178.3	16.0	0.6	179 P, 194 P, 162 P
		5000 µg	112.0	7.0	0.4	104 P, 117 P, 115 P
	<b>DMSO</b>		296.0	14.7		279, 304, 305
	<b>Untreated</b>		327.7	20.5		348, 328, 307
<b>WP2 uvrA pKM101</b>	<b>Profenofos EC (A8591B)</b>	3 µg	307.3	7.1	1.1	301, 306, 315
		10 µg	305.0	24.1	1.0	280, 307, 328
		33 µg	297.7	8.5	1.0	301, 288, 304
		100 µg	305.0	11.0	1.0	305, 316, 294
		333 µg	272.0	12.3	0.9	286, 263, 267
		1000 µg	232.0	21.1	0.8	254, 212, 230
		2500 µg	200.3	6.1	0.7	207 P, 195 P, 199 P
		5000 µg	180.7	8.1	0.6	188 P M, 182 P M, 172 P M
	<b>DMSO</b>		291.3	12.7		283, 306, 285
	<b>Untreated</b>		322.0	16.1		335, 327, 304
<b>TA 1535</b>	<b>NaN3</b>	10 µg	1040.0	73.7	80.0	1055, 960, 1105
<b>TA 1537</b>	<b>4-NOPD</b>	50 µg	78.3	4.2	5.5	83, 77, 75
<b>TA 98</b>	<b>4-NOPD</b>	10 µg	721.3	70.4	33.3	640, 761, 763
<b>TA 100</b>	<b>NaN3</b>	10 µg	1561.7	96.2	15.2	1642, 1588, 1455
<b>WP2 pKM101</b>	<b>MMS</b>	2.0 µL	2749.0	102.2	9.3	2855, 2651, 2741
<b>WP2 uvrA pKM101</b>	<b>MMS</b>	2.0 µL	2963.7	226.0	10.2	2969, 2735, 3187
Key to Positive Controls						Key to Plate Postfix Codes
NaN3	sodium azide					P Precipitate
4-NOPD	4-nitro-o-phenylene-diamine					M Manual count
MMS	methyl methane sulfonate					

**TABLE 4                      Pre-Experiment/Experiment I: 2180100 VV Plate Incorporation With Metabolic Activation**

Study Name: 2180100  
Experiment: 2180100 VV Plate  
Assay Conditions:

Study Code: ICCR 2180100  
Date Plated: 11.08.2021  
Date Counted: 18.08.2021

With metabolic activation						
Strain	Compound	Concen- tration per plate	Mean revertants per plate	Standard Deviation	Ratio treated / solvent	Individual revertant colony counts
TA 1535	Profenofos EC (A8591B)	3 µg	14.0	2.6	1.3	11, 16, 15
		10 µg	14.0	4.4	1.3	17, 16, 9
		33 µg	12.7	4.0	1.2	17, 12, 9
		100 µg	13.3	1.2	1.2	14, 12, 14
		333 µg	12.3	2.9	1.1	9, 14, 14
		1000 µg	10.7	3.5	1.0	11, 7, 14
		2500 µg	8.7	2.9	0.8	12 P, 7 P, 7 P
		5000 µg	9.7	1.2	0.9	9 P M, 11 P M, 9 P M
	DMSO		11.0	2.6		10, 14, 9
	Untreated		12.3	2.5		10, 15, 12
TA 1537	Profenofos EC (A8591B)	3 µg	18.3	3.1	1.0	19, 21, 15
		10 µg	16.0	5.0	0.9	11, 21, 16
		33 µg	16.7	2.9	0.9	20, 15, 15
		100 µg	17.0	2.0	0.9	17, 19, 15
		333 µg	11.7	2.5	0.6	14, 9, 12
		1000 µg	11.0	1.7	0.6	9, 12, 12
		2500 µg	11.3	2.1	0.6	9 P M, 13 P M, 12 P M
		5000 µg	11.0	3.0	0.6	8 P M, 11 P M, 14 P M
	DMSO		18.7	1.5		20, 17, 19
	Untreated		18.3	3.1		21, 15, 19
TA 98	Profenofos EC (A8591B)	3 µg	52.7	8.5	1.2	56, 43, 59
		10 µg	52.7	3.5	1.2	49, 53, 56
		33 µg	57.0	13.0	1.3	44, 70, 57
		100 µg	50.0	1.7	1.2	49, 52, 49
		333 µg	33.3	5.5	0.8	27, 36, 37
		1000 µg	33.3	11.0	0.8	46, 26, 28
		2500 µg	25.3	8.1	0.6	16 P, 30 P, 30 P
		5000 µg	26.3	3.1	0.6	23 P M, 29 P M, 27 P M
	DMSO		43.3	13.3		58, 40, 32
	Untreated		48.7	2.5		51, 49, 46
Key to Plate Postfix Codes						
						P      Precipitate
						M      Manual count

Study Name: 2180100  
 Experiment: 2180100 VV Plate  
 Assay Conditions:

Study Code: ICCR 2180100  
 Date Plated: 11.08.2021  
 Date Counted: 18.08.2021

**With metabolic activation**

Strain	Compound	Concentration per plate	Mean revertants per plate	Standard Deviation	Ratio treated / solvent	Individual revertant colony counts
<b>TA 100</b>	<b>Profenofos EC (A8591B)</b>	3 µg	110.0	8.7	1.1	120, 105, 105
		10 µg	109.3	11.0	1.1	122, 104, 102
		33 µg	119.7	3.1	1.2	119, 117, 123
		100 µg	109.3	5.5	1.1	109, 115, 104
		333 µg	97.0	5.6	1.0	98, 91, 102
		1000 µg	76.3	12.7	0.8	69, 91, 69
		2500 µg	39.7	10.4	0.4	43 P, 28 P, 48 P
		5000 µg	29.3	7.4	0.3	32 P M, 21 P M, 35 P M
	<b>DMSO</b>		100.7	19.6		106, 79, 117
	<b>Untreated</b>		116.7	13.3		102, 128, 120
<b>WP2 pKM101</b>	<b>Profenofos EC (A8591B)</b>	3 µg	326.3	19.1	1.1	332, 342, 305
		10 µg	300.3	6.1	1.0	295, 299, 307
		33 µg	278.7	24.7	0.9	307, 267, 262
		100 µg	297.3	9.8	1.0	286, 303, 303
		333 µg	254.0	28.2	0.8	275, 265, 222
		1000 µg	187.3	15.3	0.6	205, 178, 179
		2500 µg	160.3	21.1	0.5	163 P, 138 P, 180 P
		5000 µg	111.7	5.0	0.4	107 P M, 117 P M, 111 P M
	<b>DMSO</b>		310.3	34.2		274, 342, 315
	<b>Untreated</b>		340.3	5.9		338, 347, 336
<b>WP2 uvrA pKM101</b>	<b>Profenofos EC (A8591B)</b>	3 µg	362.7	6.7	1.0	370, 357, 361
		10 µg	363.3	9.1	1.0	370, 367, 353
		33 µg	343.0	21.2	0.9	359, 351, 319
		100 µg	309.3	30.0	0.8	303, 342, 283
		333 µg	339.3	24.8	0.9	335, 366, 317
		1000 µg	276.0	36.3	0.8	245, 316, 267
		2500 µg	234.7	7.6	0.6	226 P, 238 P, 240 P
		5000 µg	162.7	12.9	0.4	148 P M, 172 P M, 168 P M
	<b>DMSO</b>		365.7	3.2		362, 368, 367
	<b>Untreated</b>		361.3	24.6		359, 387, 338
<b>TA 1535</b>	<b>2-AA</b>	2.5 µg	296.7	8.0	27.0	296, 305, 289
<b>TA 1537</b>	<b>2-AA</b>	2.5 µg	385.3	48.3	20.6	438, 375, 343
<b>TA 98</b>	<b>2-AA</b>	2.5 µg	2350.3	78.1	54.2	2439, 2292, 2320
<b>TA 100</b>	<b>2-AA</b>	2.5 µg	3439.7	40.2	34.2	3414, 3486, 3419
<b>WP2 pKM101</b>	<b>2-AA</b>	10.0 µg	952.3	69.7	3.1	958, 880, 1019
<b>WP2 uvrA pKM101</b>	<b>2-AA</b>	10.0 µg	1670.0	45.2	4.6	1635, 1721, 1654

Key to Positive Controls

2-AA 2-aminoanthracene

Key to Plate Postfix Codes

P Precipitate  
 M Manual count

**TABLE 5 Experiment II: 2180100 HV2 Pre Incubation Without Metabolic Activation**

Study Name: 2180100  
Experiment: 2180100 HV2 Pre  
Assay Conditions:

Study Code: ICCR 2180100  
Date Plated: 01.09.2021  
Date Counted: 08.09.2021

Without metabolic activation						
Strain	Compound	Concentration per plate	Mean revertants per plate	Standard Deviation	Ratio treated / solvent	Individual revertant colony counts
TA 1535	Profenofos EC (A8591B)	10 µg	13.0	3.5	0.9	9, 15, 15
		33 µg	11.3	0.6	0.8	12, 11, 11
		100 µg	14.3	4.5	1.0	19, 14, 10
		333 µg	12.0	3.5	0.9	16, 10, 10
		1000 µg	8.3	2.1	0.6	6, 9, 10
		2500 µg	8.7	1.5	0.6	10 P, 7 P, 9 P
		5000 µg	9.0	0.0	0.6	9 P M, 9 P M, 9 P M
	DMSO		14.0	2.0		12, 14, 16
	Untreated		15.7	4.7		14, 21, 12
TA 1537	Profenofos EC (A8591B)	10 µg	10.3	3.5	0.9	7, 10, 14
		33 µg	10.3	1.5	0.9	12, 10, 9
		100 µg	11.7	2.5	1.0	14, 12, 9
		333 µg	11.0	2.6	0.9	10, 14, 9
		1000 µg	8.3	2.9	0.7	5, 10, 10
		2500 µg	8.3	1.2	0.7	9 P, 9 P, 7 P
		5000 µg	5.7	2.1	0.5	4 P M, 5 P M, 8 P M
	DMSO		11.7	2.5		14, 12, 9
	Untreated		9.7	2.3		11, 11, 7
TA 98	Profenofos EC (A8591B)	10 µg	25.7	4.5	1.1	21, 30, 26
		33 µg	28.3	8.7	1.2	26, 21, 38
		100 µg	21.0	6.2	0.9	23, 26, 14
		333 µg	19.7	4.6	0.9	17, 25, 17
		1000 µg	22.7	2.1	1.0	21, 25, 22
		2500 µg	20.0	0.0	0.9	20 P, 20 P, 20 P
		5000 µg	12.7	2.5	0.6	15 P M, 13 P M, 10 P M
	DMSO		22.7	2.1		22, 25, 21
	Untreated		32.7	5.9		35, 26, 37
TA 100	Profenofos EC (A8591B)	10 µg	103.3	9.0	1.0	94, 104, 112
		33 µg	93.0	7.0	0.9	101, 88, 90
		100 µg	81.7	11.0	0.8	94, 78, 73
		333 µg	83.7	9.0	0.8	79, 78, 94
		1000 µg	76.0	7.9	0.7	79, 82, 67
		2500 µg	38.0	8.7	0.4	28 P, 43 P, 43 P
		5000 µg	18.0	3.0	0.2	21 P M, 18 P M, 15 P M
	DMSO		104.7	6.4		100, 112, 102
	Untreated		93.7	12.0		93, 106, 82

Key to Plate Postfix Codes

P Precipitate  
M Manual count

Study Name: 2180100  
 Experiment: 2180100 HV2 Pre  
 Assay Conditions:

Study Code: ICCR 2180100  
 Date Plated: 01.09.2021  
 Date Counted: 08.09.2021

**Without metabolic activation**

Strain	Compound	Concentration per plate	Mean revertants per plate	Standard Deviation	Ratio treated / solvent	Individual revertant colony counts
<b>WP2 pKM101</b>	<b>Profenofos EC (A8591B)</b>	10 µg	259.7	5.1	1.0	264, 261, 254
		33 µg	263.7	11.0	1.0	263, 253, 275
		100 µg	177.7	22.8	0.7	165, 204, 164
		333 µg	155.3	10.2	0.6	148, 151, 167
		1000 µg	119.0	21.7	0.5	116, 99, 142
		2500 µg	95.7	17.9	0.4	75 P, 106 P, 106 P
		5000 µg	64.3	7.8	0.2	73 P, 62 P, 58 P
	<b>DMSO</b>		262.7	30.6		274, 286, 228
	<b>Untreated</b>		324.0	25.2		351, 320, 301
<b>WP2 uvrA pKM101</b>	<b>Profenofos EC (A8591B)</b>	10 µg	268.3	21.1	1.0	286, 274, 245
		33 µg	258.0	7.9	0.9	249, 264, 261
		100 µg	249.7	3.1	0.9	249, 247, 253
		333 µg	256.0	13.9	0.9	263, 240, 265
		1000 µg	181.0	12.1	0.7	179, 194, 170
		2500 µg	178.7	15.0	0.7	178 P, 164 P, 194 P
		5000 µg	132.7	2.5	0.5	130 P, 133 P, 135 P
	<b>DMSO</b>		274.7	5.5		275, 269, 280
	<b>Untreated</b>		313.7	10.6		325, 304, 312
<b>TA 1535</b>	<b>NaN3</b>	10 µg	938.0	44.2	67.0	987, 901, 926
<b>TA 1537</b>	<b>4-NOPD</b>	50 µg	100.0	14.2	8.6	116, 95, 89
<b>TA 98</b>	<b>4-NOPD</b>	10 µg	777.0	41.1	34.3	820, 773, 738
<b>TA 100</b>	<b>NaN3</b>	10 µg	1425.3	142.9	13.6	1574, 1289, 1413
<b>WP2 pKM101</b>	<b>MMS</b>	2.0 µL	2538.3	35.2	9.7	2499, 2549, 2567
<b>WP2 uvrA pKM101</b>	<b>MMS</b>	2.0 µL	2205.3	130.8	8.0	2154, 2354, 2108
Key to Positive Controls						Key to Plate Postfix Codes
NaN3	sodium azide					P Precipitate
4-NOPD	4-nitro-o-phenylene-diamine					M Manual count
MMS	methyl methane sulfonate					



**TABLE 6 Experiment II: 2180100 HV2 Pre Incubation With Metabolic Activation**

Study Name: 2180100  
Experiment: 2180100 HV2 Pre  
Assay Conditions:

Study Code: ICCR 2180100  
Date Plated: 01.09.2021  
Date Counted: 08.09.2021

**With metabolic activation**

Strain	Compound	Concen- tration per plate	Mean revertants per plate	Standard Deviation	Ratio treated / solvent	Individual revertant colony counts
TA 1535	Profenofos EC (A8591B)	10 µg	14.0	2.0	1.2	14, 16, 12
		33 µg	13.7	2.5	1.2	14, 16, 11
		100 µg	10.7	0.6	0.9	10, 11, 11
		333 µg	12.0	2.6	1.1	10, 15, 11
		1000 µg	13.7	3.2	1.2	16, 10, 15
		2500 µg	14.0	3.5	1.2	16 P, 16 P, 10 P
		5000 µg	9.0	3.0	0.8	12 P M, 9 P M, 6 P M
	DMSO		11.3	2.3		10, 10, 14
	Untreated		10.7	1.2		10, 12, 10
TA 1537	Profenofos EC (A8591B)	10 µg	13.7	3.2	0.8	10, 15, 16
		33 µg	13.3	3.2	0.8	11, 12, 17
		100 µg	13.0	1.7	0.8	11, 14, 14
		333 µg	12.3	1.5	0.8	11, 14, 12
		1000 µg	13.7	2.3	0.8	15, 15, 11
		2500 µg	13.0	1.7	0.8	12 P, 12 P, 15 P
		5000 µg	9.7	1.2	0.6	11 P M, 9 P M, 9 P M
	DMSO		16.3	0.6		16, 16, 17
	Untreated		13.0	2.6		11, 16, 12
TA 98	Profenofos EC (A8591B)	10 µg	44.3	12.7	1.0	30, 49, 54
		33 µg	45.7	4.2	1.0	49, 47, 41
		100 µg	48.3	5.0	1.1	43, 53, 49
		333 µg	35.3	7.4	0.8	38, 27, 41
		1000 µg	28.3	4.7	0.6	32, 30, 23
		2500 µg	28.0	2.0	0.6	26 P, 28 P, 30 P
		5000 µg	20.7	6.0	0.5	27 P M, 20 P M, 15 P M
	DMSO		45.0	9.6		56, 41, 38
	Untreated		39.0	9.5		38, 49, 30
TA 100	Profenofos EC (A8591B)	10 µg	104.0	16.0	1.0	88, 104, 120
		33 µg	113.7	15.0	1.1	131, 106, 104
		100 µg	103.3	16.6	1.0	121, 101, 88
		333 µg	86.3	6.4	0.9	79, 89, 91
		1000 µg	87.0	8.5	0.9	95, 88, 78
		2500 µg	39.0	6.9	0.4	35 P, 47 P, 35 P
		5000 µg	25.0	3.6	0.2	28 P M, 21 P M, 26 P M
	DMSO		101.0	9.8		90, 104, 109
	Untreated		98.7	16.8		84, 95, 117

Key to Plate Postfix Codes

P Precipitate  
M Manual count

Study Name: 2180100  
 Experiment: 2180100 HV2 Pre  
 Assay Conditions:

Study Code: ICCR 2180100  
 Date Plated: 01.09.2021  
 Date Counted: 08.09.2021

**With metabolic activation**

Strain	Compound	Concentration per plate	Mean revertants per plate	Standard Deviation	Ratio treated / solvent	Individual revertant colony counts
<b>WP2 pKM101</b>	<b>Profenofos EC (A8591B)</b>	10 µg	291.0	13.9	0.9	298, 275, 300
		33 µg	272.3	30.2	0.9	307, 252, 258
		100 µg	283.7	29.6	0.9	258, 316, 277
		333 µg	219.3	27.2	0.7	237, 233, 188
		1000 µg	170.0	10.4	0.5	165, 182, 163
		2500 µg	174.7	27.8	0.6	143 P, 186 P, 195 P
		5000 µg	118.7	7.6	0.4	112 P, 127 P, 117 P
	<b>DMSO</b>		311.7	26.6		284, 337, 314
	<b>Untreated</b>		319.3	21.8		300, 315, 343
<b>WP2 uvrA pKM101</b>	<b>Profenofos EC (A8591B)</b>	10 µg	368.7	27.0	1.0	368, 342, 396
		33 µg	360.0	13.7	1.0	372, 363, 345
		100 µg	355.7	22.4	1.0	380, 336, 351
		333 µg	358.3	25.5	1.0	384, 333, 358
		1000 µg	269.7	18.4	0.7	290, 265, 254
		2500 µg	221.3	7.8	0.6	230 P, 215 P, 219 P
		5000 µg	164.0	6.0	0.4	170 P M, 158 P M, 164 P M
	<b>DMSO</b>		372.7	15.6		387, 356, 375
	<b>Untreated</b>		392.3	8.3		399, 395, 383
<b>TA 1535</b>	<b>2-AA</b>	2.5 µg	312.7	23.2	27.6	324, 286, 328
<b>TA 1537</b>	<b>2-AA</b>	2.5 µg	307.0	43.1	18.8	348, 311, 262
<b>TA 98</b>	<b>2-AA</b>	2.5 µg	2288.3	267.8	50.9	2082, 2591, 2192
<b>TA 100</b>	<b>2-AA</b>	2.5 µg	2502.0	175.3	24.8	2696, 2355, 2455
<b>WP2 pKM101</b>	<b>2-AA</b>	10.0 µg	959.3	34.1	3.1	992, 962, 924
<b>WP2 uvrA pKM101</b>	<b>2-AA</b>	10.0 µg	1451.7	89.9	3.9	1509, 1498, 1348
Key to Positive Controls						Key to Plate Postfix Codes
2-AA	2-aminoanthracene					P Precipitate M Manual count

## **APPENDICES SECTION**

## APPENDIX 1 Historical Control Data

These data represent the laboratory's historical control data from July 2018 until July 2020 representing approx. 600 experiments (WP2 pKM101, WP2 uvrA pKM101 the historical data are based on approx. 80 experiments).

The positive controls that used to compile the historical positive control data correspond to the positive control substances described in Methods; section 3.2.2 (Positive control substances).

Strain		without S9 mix				with S9 mix			
		Mean	SD	Min	Max	Mean	SD	Min	Max
TA 1535	Solvent control	12	2.6	7	22	13	2.5	7	24
	Untreated control	12	2.9	6	26	13	2.8	7	23
	Positive control	1116	141.3	340	1612	346	72.1	170	736
TA1537	Solvent control	11	2.4	6	20	14	2.8	7	28
	Untreated control	11	2.8	5	22	14	3.2	7	30
	Positive control	83	22.1	48	400	286	98.7	82	630
TA 98	Solvent control	28	4.9	13	46	38	6.4	12	62
	Untreated control	29	5.0	14	48	41	6.8	14	64
	Positive control	421	91.2	216	1218	3275	774.9	322	5699
TA 100	Solvent control	127	30.7	63	214	131	30.0	72	214
	Untreated control	135	35.7	64	233	140	34.4	68	217
	Positive control	1759	273.4	511	2588	3566	837.6	553	5444
WP2 pKM 101	Solvent control	248	31.7	171	299	266	33.0	205	315
	Untreated control	269	26.6	212	346	299	28.2	233	345
	Positive control	3343	428.4	2332	4653	1092	257.8	933	2781
WP2uvrA pKM 101	Solvent control	322	31.6	248	388	375	38.5	287	466
	Untreated control	346	28.2	279	403	393	32.6	313	480
	Positive control	3176	468.5	2021	4717	1897	183.2	1270	2464

Mean = mean value of revertants/plate

SD = standard deviation

Min = minimal value

Max = maximal value

## APPENDIX 2 Copy of GLP Certificate



### Gute Laborpraxis/Good Laboratory Practice

### GLP-Bescheinigung/Statement of GLP Compliance

(gemäß/according to § 19b Abs. 1 Chemikaliengesetz)



Eine GLP-Inspektion zur Überwachung der Einhaltung der GLP-Grundsätze gemäß Chemikaliengesetz bzw. Richtlinie 2004/9/EG wurde durchgeführt in

Assessment of conformity with GLP according to Chemikaliengesetz and Directive 2004/9/EEC at:

☒ Prüfeinrichtung/Test facility ☐ Prüfstandort/Test site

**ICCR-Roßdorf GmbH**  
Institute for Competent Contract Research  
In den Leppsteinswiesen 19  
64380 Roßdorf

(Unverwechselbare Bezeichnung und Adresse/Unequivocal name and address)

### Prüfungen nach Kategorien/Areas of Expertise

(gemäß/according ChemVwV-GLP Nr. 5.3/OECD guidance)

**2** Prüfungen zur Bestimmung der toxikologischen Eigenschaften

**2** Toxicity studies

**3** Prüfungen zur Bestimmung der erbgutverändernden Eigenschaften (in vitro und in vivo)

**3** Mutagenicity studies

**8** Analytische Prüfungen an biologischen Materialien

**8** Analytical and clinical chemistry testing

**22.11.2018, 21.02.2019, 12. bis 14.03.2019**  
Datum der Inspektion/Date of Inspection  
(Tag Monat Jahr/day month year)

Die genannte Prüfeinrichtung befindet sich im nationalen GLP-Überwachungsverfahren und wird regelmäßig auf Einhaltung der GLP-Grundsätze überwacht.

The above mentioned test facility is included in the national GLP Compliance Programme and is inspected on a regular basis.

Auf der Grundlage des Inspektionsberichtes wird hiermit bestätigt, dass in dieser Prüfeinrichtung die oben genannten Prüfungen unter Einhaltung der GLP- Grundsätze durchgeführt werden können.

Based on the inspection report it can be confirmed, that this test facility is able to conduct the aforementioned studies in compliance with the Principles of GLP.

Im Auftrag

Dr. Astrid Brandt, Referentin, Wiesbaden, den **23. Oktober 2019**  
(Name und Funktion der verantwortlichen Person/  
Name and function of responsible person)



Hessisches Ministerium für Umwelt, Klimaschutz, Landwirtschaft und Verbraucherschutz,  
Mainzer Straße 80, D 65189 Wiesbaden  
(Name und Adresse der GLP-Überwachungsbehörde/Name and address of the GLP Monitoring Authority)

English name and address of the GLP Monitoring Authority: Hessian Ministry for Environment, Climate Protection, Agriculture and Consumer Protection; Department II 10; P.O. Box 31 09; 65189 Wiesbaden

Translation of seal inscription: Hessian Ministry for Environment, Climate Protection, Agriculture and Consumer Protection

## APPENDIX 3 Certificate of S9



### CERTIFICATE

**ICCR-Roßdorf S9 Preparation Lot No. 200521**

**Date of preparation: May 20, 2021**

**Release date: June 07, 2021**

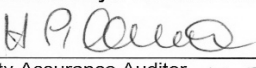
**Protein assay: 30.2 mg protein / ml S9**

**Sterility: 0 colonies / ml S9 on glucose-minimal-agar**

**Salmonella typhimurium assay (AMES-test)**

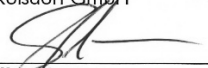
Treatment	µl S9 / plate	number of revertants in TA 98
negative	0	29
control	100	68
10 µg/plate	0	130
2-Aminoanthracene	150	2344
10 µg/plate	0	34
Benzo(a)pyrene	100	105

The S9 was obtained from the livers of male Wistar rats which received triple treatments of 80 mg / kg body weight Phenobarbital and  $\beta$ -Naphthoflavone orally on consecutive days. The livers were prepared 24 hours after the last treatment.

  
Quality Assurance Auditor  
ICCR-Roßdorf GmbH **H. Pilawa**

**18. JUNI 2021**

Date

  
Dr. Steffen Naumann  
Study Director  
ICCR-Roßdorf GmbH

**18. JUNI 2021**

Date

**ICCR-Roßdorf GmbH**  
In den Leppsteinswiesen 19, 64380 Roßdorf, Deutschland  
T +49 6154 8070 F +49 6154 83399  
Registergericht Darmstadt, HRB 6837, USt-ID DE812339696  
Geschäftsführer: Dr. Markus Schulz

*SOP Origin TS-SOP S9\_23*

## APPENDIX 4 Certificate of Analysis



ALS Laboratórios LS Ltda.  
Rua Fábria, 59 – CEP: 05051-030  
São Paulo, SP - Brazil

SYNGENTA PROTEÇÃO DE CULTIVOS Ltda.  
Rua Doutor Rubens Gomes Bueno nº 691,  
11º andar, Torre Sigma  
CEP 04730-000 – Bairro Várzea de Baixo  
São Paulo-SP – Brazil

### Certificate of Analysis

**A8591B**  
**Profenofos EC (960)**  
**RAN001-099-019**

Batch Identification	RAN001-099-019
Product Code	A8591B
Other Product Code(s)	A8591; CGA15324 EC (960); EXF23490E
EUP number	514/2020 Expiry date: 26/02/2023
Received on	12 May 2021
Source	Syngenta Proteção de Cultivos Ltda. Rodovia Professor Zeferino Vaz, SP 332, s/nº, km 127,5 – Bairro Santa Terezinha, CEP 13148-915 – Paulínia – SP – Brasil

#### Chemical Analysis (Active Ingredients Content)

– Content of Profenofos *	73.12 % w/w corresponding to 969.03 g/L
---------------------------	---

The Active Ingredient content is within the FAO limits.

Methodology used for Characterization: CG-FID (SF-1135/1)

#### Physical Analysis

– Density *	1.3253 g/cm <sup>3</sup>
-------------	--------------------------

#### Stability:

– Storage Temperature	<30°C
– Recertification Date	06 May 2023

If stored under the conditions given above, this test item can be considered stable until the recertification date is reached.

This Certificate of Analysis summarizes data which originates either from a single study or from several individual studies. Tests marked with an asterisk (\*) have been conducted in compliance with GLP. All original raw data, including any storage medium for electronically recorded data, documentation, the signed study plan, the protocol amendments, the final report and a sample of the test item will be retained in the GLP Archives at ALS Laboratórios LS Ltda.

Study number of batch characterization: 25926/2021CC

Authorization: 26 May 2021

*Victor F. G. da Silva*  
Victor Ferreira Gomes da Silva  
ALS Laboratórios LS Ltda.